



DEVELOPMENT AND EVALUATION OF A MULTIPLEX CONVENTIONAL RT-PCR FOR SIMULTANEOUS DETECTION AND TYPING OF FMDV IN WEST AFRICA

Kamila Gorna, Gaye Laye Diop, Marianne Diop, Eve Laloy, Modou
Moustapha Lo, Anthony Relmy, Aurore Romey, Momar Talla Seck, Hussaini
Ularamu, Stéphan Zientara, et al.

► To cite this version:

Kamila Gorna, Gaye Laye Diop, Marianne Diop, Eve Laloy, Modou Moustapha Lo, et al.. DEVELOPMENT AND EVALUATION OF A MULTIPLEX CONVENTIONAL RT-PCR FOR SIMULTANEOUS DETECTION AND TYPING OF FMDV IN WEST AFRICA. EuFMD Open Session 2018, Oct 2018, Borgo Egnazzia, Italy. anses-04441092

HAL Id: anses-04441092

<https://anses.hal.science/anses-04441092>

Submitted on 6 Feb 2024

HAL is a multi-disciplinary open access archive for the deposit and dissemination of scientific research documents, whether they are published or not. The documents may come from teaching and research institutions in France or abroad, or from public or private research centers.

L'archive ouverte pluridisciplinaire **HAL**, est destinée au dépôt et à la diffusion de documents scientifiques de niveau recherche, publiés ou non, émanant des établissements d'enseignement et de recherche français ou étrangers, des laboratoires publics ou privés.

CONTEXT

West Africa : circulation of O/A/SAT1/SAT2 types ► Pool 5




- ✓ FMD is endemic in West Africa:
- Continuous FMDV circulation
- ✓ Lack of diagnostic capacities & under-reporting:
- Lack of information on FMD dynamics
- ✓ Difficulties for implementation of control strategies
- risk analysis of FMDV incursion/spread?
- ✓ Towards the goal of FMD eradication...
- Knowledge definitely required on FMDV strains circulating in West Africa




AIM

➔ To develop and evaluate a conventional multiplex RT-PCR for both detection and typing of FMD viruses circulating in West Africa


Target sequences

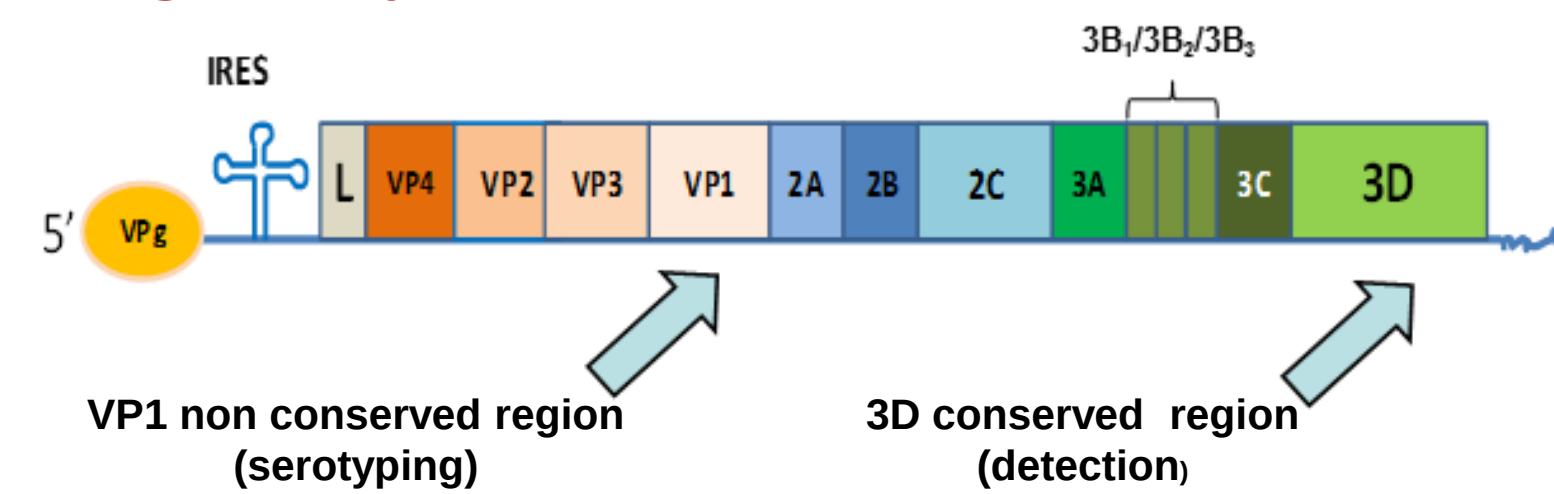

Genome analysis & Design of primers

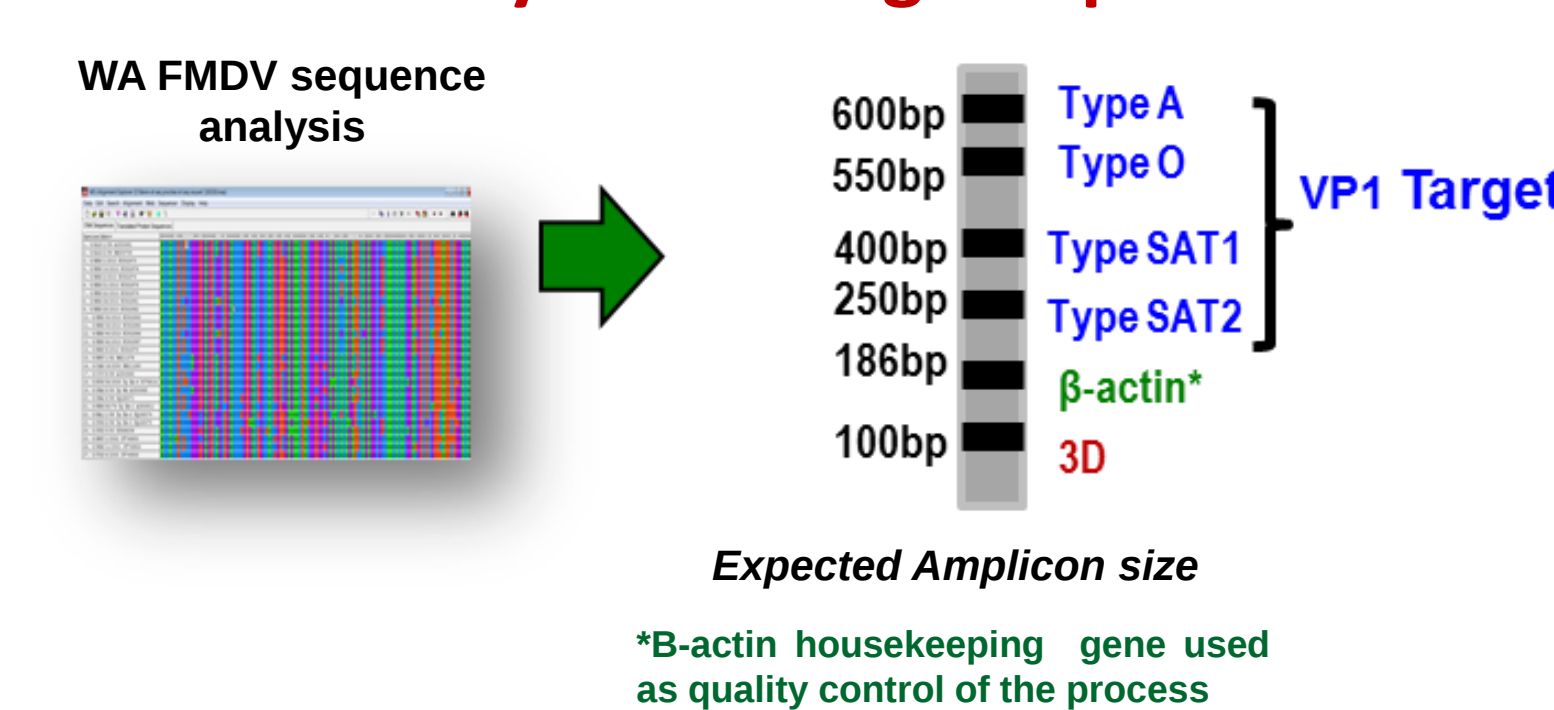


Evaluation Steps

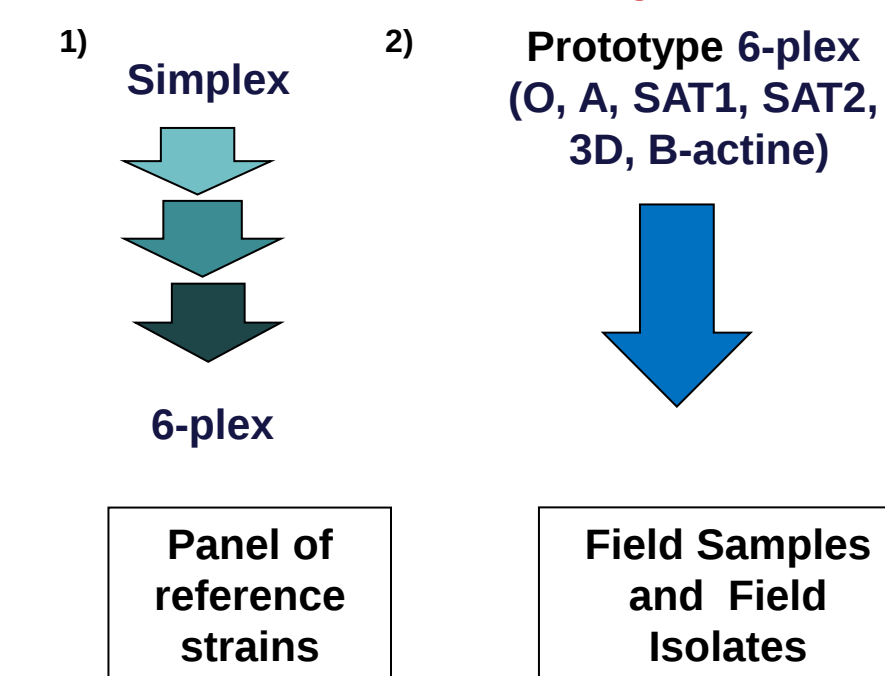

Target specificity & amplicon sizes allowing easy visual interpretation of RT-PCR profiles using agarose gel electrophoresis

6 PLEX RT-PCR PROTOTYPE (O/A/SAT1/SAT2/3D/B-ACTIN)

► Developed using reference strains and field isolates

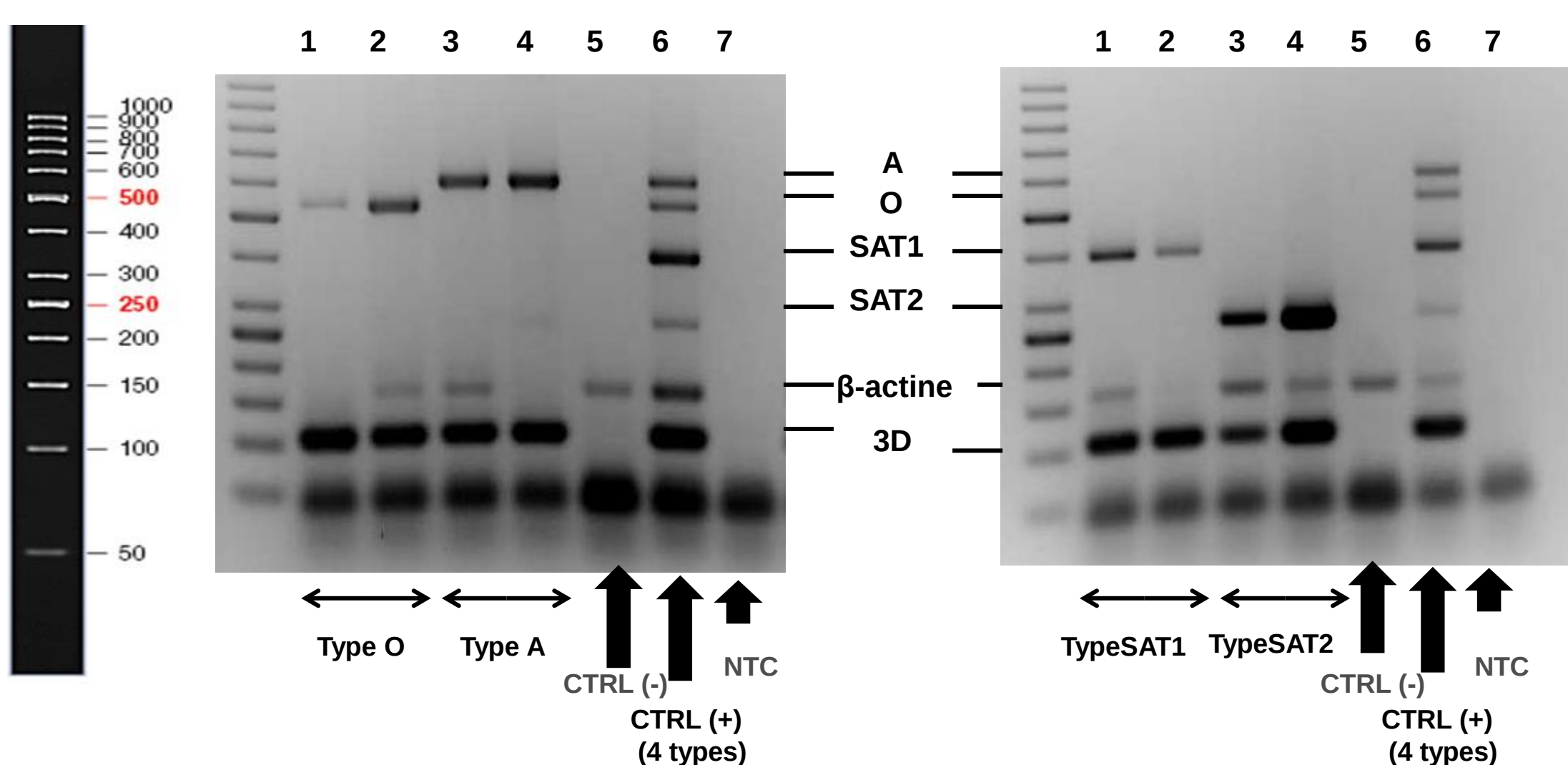

Analytical Specificity of RT-PCR 6-plex

Strain serotype	n	3D	Typing
Type O	7	7/7	7/7
Originated from West Africa	2	2/2	2/2
Other isolates	5	5/5	5/5
Type A	8	8/8	8/8
Originated from West Africa	2	2/2	2/2
Other isolates	6	6/6	6/6
Type SAT1 strains (Kenilbot)	2	2/2	2/2
Type SAT2 strains (Libitigayim)	3	3/3	3/3
Type SAT3 strains (Zim)	1	1/1	0/1
Type C	1	1/1	0/1
Type Asia-1	2	2/2	0/2
EVDF*	3	0/3	0/3

* Swine vesicular disease virus (i.e. other picornavirus) (differential diagnostic of FMD)

► Diagnostic sensitivity & specificity

Diagnostic sensitivity (on positive samples type O & A)


Positive field samples and the corresponding isolates (West Africa; Benin 2010)

Sample nature	Type	Real-time RT-PCR (2D Target)	RT-PCR 6-plex		RT-PCR 3-plex		Simplex (VP1 typing O or A)
			Detection (2D)	Typing (VP1 O or A)	Detection (2D)	Typing (VP1 O or A)	
Field Samples	O	28/28	28/28	8/28	28/28	13/28	28/28
Field Isolates	O	28/28	28/28	24/28	28/28	28/28	28/28
Field Samples	A	9/9	9/9	6/9	9/9	6/9	9/9
Field Isolates	A	9/9	9/9	9/9	9/9	9/9	9/9

Diagnostic specificity (negative samples)

- ✓ Negative Fields samples (n=24) /Bovine, Swine, Ovine, Caprine/ (France 2013-14)
- ✓ None amplification was observed for targets 3D and VP1.
- ✓ B-actin amplification was observed for all samples

EVALUATION OF THIS PROTOTYPE ON FIELD SAMPLES THROUGH COLLABORATIONS IN ENDEMIC COUNTRIES

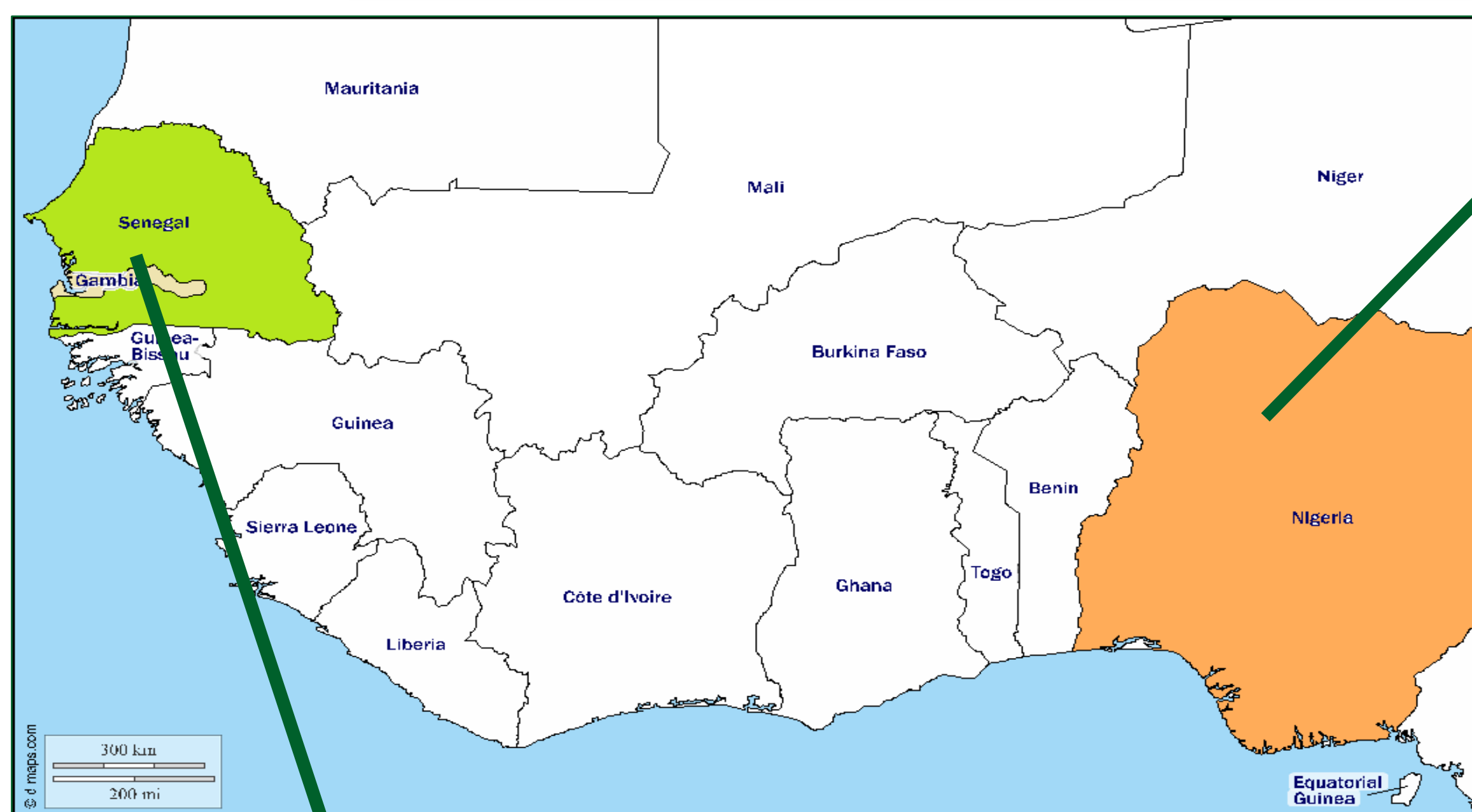

► Evaluation on field samples (Senegal, ISRA/LNERV)

- 51 samples
- Collected between 2008-2014
- From Senegal & Gambia


+ wellcome of a trainee (2 weeks) from LNERV

► Evaluation on field samples & corresponding isolates (Nigeria, NVRI)

- 40 epithelium samples (38 bovine/2 ovine)
- Collected between 2011-2016
- From 6 states

S/No	Sample ID	RT-PCR	VI	Serotype	OD	1st res.	2nd res.
1	PL/BA/08/106	Pos	2nd P	SAT2	0.201	SAT2	/
2	PL/BA/08/108	Pos	2nd P	SAT2	0.306	SAT2	/
3	PL/DN/00/6/E	Pos	2nd P	O	1.336	O/SAT2	/
4	PL/DN/00/6/E	Pos	2nd P	O	0.346	O	/
5	PL/DN/00/6/E	Pos	2nd P	SAT2	0.339	SAT2	/
6	BAU/10/13	Pos	2nd P	A	0.776	A/O	/
7	BAU/10/13	Pos	2nd P	O	0.955	O	/
8	BAU/10/13	Pos	2nd P	A	0.786	O	/
9	BAU/10/13	Pos	2nd P	A	0.36	A	/
10	PL/KA/06/04/1	Pos	2nd P	O	1.137	O	/
11	PL/KA/06/04/18-1	Pos	2nd P	O	1.217	O	/
12	PL/KA/06/04/18-2	Pos	2nd P	O	0.954	O	/
13	KD/BAU/01/12	Pos	2nd P	O	0.868	O	/
14	BN/G/16	Pos	2nd P	O	1.281	O	/
15	PL/BA/12/16	Pos	2nd P	O	1.289	O	/
16	PL/BA/12/16	Pos	2nd P	O	1.249	O	/
17	PL/BA/12/16	Pos	2nd P	SAT2	0.523	SAT2	/
18	PL/PKN/01/16	Pos	2nd P	O	1.212	O/low	/
19	PL/PKN/01/16	Pos	2nd P	O	1.169	O	/
20	PL/JS/06/01/16	Pos	2nd P	O	0.735	O	/
21	PL/JS/06/01/16	Pos	2nd P	O	1.420	O	/
22	PL/KARA/15/04/17-1	Pos	2nd P	O	0.971	O	/
23	PL/KARA/15/04/17-1	Pos	2nd P	O	1.642	O	/
24	PL/JS/KARA/12/16/1	Pos	2nd P	O	1.025	O	/
25	PL/JS/KARA/12/16/1	Pos	2nd P	A/low	0.169	SAT2 low / O/SAT2	/
26	BA/19/1	Pos	2nd P	A/low	0.162	A	/
27	BA/19/1	Pos	2nd P	A/low	0.129	A	/
28	PL/TN/01/11	Pos	2nd P	O	0.630	A/low	/
29	PL/TN/01/11	Pos	2nd P	SAT2	1.025	SAT2	/
30	PL/78/15	Pos	2nd P	O	0.262	O	/
31	PL/80/15	Pos	2nd P	O	0.262	O	/
32	PL/32/15	Pos	2nd P	O	0.262	O	/
33	PL/57/15	Pos	2nd P	SAT1	0.226	SAT1 (low)	/
34	KW/12/14	Pos	2nd P	A	0.641	A	/
35	PL/BA/03/08	Pos	2nd P	A	0.566	A	/
36	PL/BA/03/08	Pos	2nd P	A	0.211	A	/
37	AD/GMB/03	Pos	2nd P	SAT2	0.435	SAT2	/
38	FMD/19/28	Pos	2nd P	O	0.587	O	/
39	BN/FMD/2011/06	Pos	2nd P	O	0.587	O	/
40	BN/FMD/2011/06	Pos	2nd P	O	0.490	O	/

Results (6-plex vs Ag-Elisa on isolates)

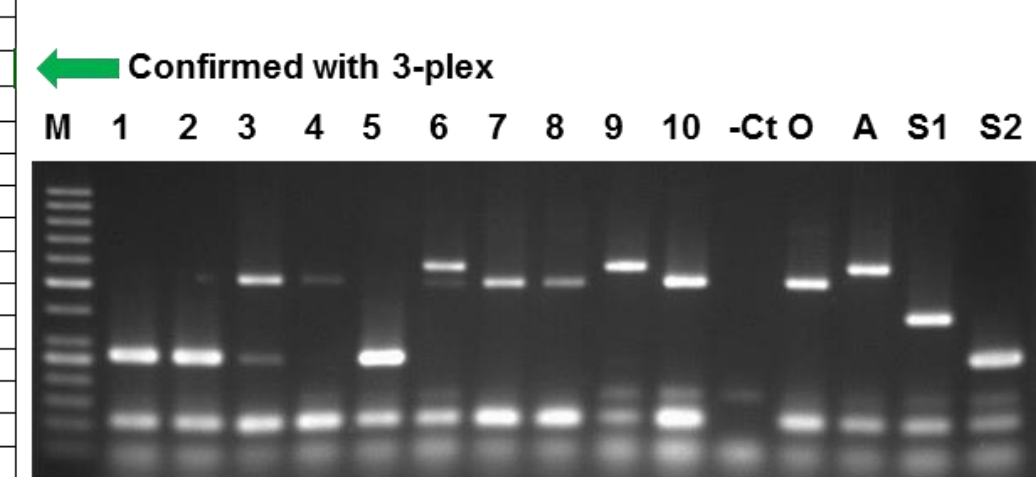

- Confirmed with 3-plex
- Sample Typed as SAT1 with 6-plex (new primers added) ✓
- Confirmed with 3-plex
- Sample Typed as SAT1 with 6-plex (new primers For/Rev added) ✓
- Samples Typed as SAT1 with 6-plex (new primers For/Rev added) ✓
- Samples Typed as A with 6-plex (new primers Rev added) ✓
- Sample Typed as SAT1 with 6-plex (new primers For/Rev added) ✓

- ✓ All isolates typed using 6plex (37/40) or 3plex confirmation (3/40)
- ✓ Last results: 22 field samples typed using 6plex (including 4 whose isolates not typed by Ag-Elisa)

- Sequencing and optimization of SAT2 / O typing
- Test of the optimized RT-PCR on the original samples and corresponding isolates
- Evaluation of the RT-PCR at NVRI and LNERV/ISRA on field samples (in progress)

CONCLUSION & PERSPECTIVES

A PROMISING CONVENTIONAL MULTIPLEX RT-PCR METHOD

- Allowing simultaneous detection and typing of FMD viruses
- ✓ Evaluated on reference strains, WA field samples & corresponding isolates
- ✓ Specific
- ✓ Flexible (triplex, 4-plex, 5-plex, 6-plex as required)
- Easy to implement in diagnostic laboratories
- ✓ No need of expensive laboratory equipment and reagents
- ✓ Could be used for early and rapid detection
- ✓ Currently used in Nigeria and Senegal (technological transfer and primer mixes provided)
- ✓ Ready-to-used mastermix under study

Results of the RTqPCR 6-plex assay (on the 27 samples 3D positive by RTqPCR)

Sample	Well number	Pan FA (3D)	B-actin	Typing
Ladder	1	/	/	/
SEN 1	2	POS	NEG	NEG
SEN 2	3	POS	NEG	NEG
SEN 3	4	POS	NEG	NEG
SEN 4	5	POS	POS	SAT2 (low)
SEN 5	6	POS	NEG	NEG
SEN 6	7	POS	NEG	SAT2
SEN 7	8	POS	POS	O
SEN 7B	9	POS	NEG	SAT2 (low)
TG/O/SAT1/SAT2	10	POS	POS (low)	POS

sample type	analysed	RtqPCR Duplex pan-FA		Rt-PCR 6 plex	
		3D positive	beta actin	3D positive	Typing
Epithelium sample (bovine)	15	14	15	14 (Ct<25)	9 (Ct<27)
Epithelium sample (swine)	2	1	2	1 (Ct<26)	0
Tongue swab (bovine)	14	6	10	1 (Ct<22)	3 (Ct<30)
Vesicles swab (bovine)	6	4	4	2 (Ct<24)	2 (Ct<27)
Mouth swab (bovine)	2	2	2	2 (Ct<15)	1 (Ct<27)
Nasal swab (goat)	7	0	7	0	5 (Ct<27)
Blood sample (bovine)	5	0	5	ND	ND
Total	51	27	45	20	7 out of 20

- ✓ 20/27 samples detected positive and 7/20 typed using 6-plex
- ✓ Overall, 10/20 samples typed (6plex + 3plex confirmation)